


CBC-ST Plus

HEMATOLOGY CONTROLS

CONTROL

ASSAY VALUES AND EXPECTED RANGES

QCP DATA MONTHS: **JUNE, JULY, AUGUST**

LOT	ST0623
	2023-09-05

Instrument	Parameter	Low	Normal	High
		LOT ST0623L	LOT ST0623N	LOT ST0623H
SYSMEX™ KX-21/KX-21N	WBC × 10 ³ /μL	2.7 ± 0.4	7.9 ± 0.8	19.6 ± 2.4
	RBC × 10 ⁶ /μL	2.53 ± 0.16	4.14 ± 0.22	5.85 ± 0.28
	HGB g/dL	6.3 ± 0.4	12.0 ± 0.5	18.6 ± 0.7
	HGB g/L	63 ± 4	120 ± 5	186 ± 7
	HGB mmol/L	3.9 ± 0.2	7.4 ± 0.3	11.5 ± 0.5
	HCT %	18.2 ± 1.8	33.3 ± 2.4	52.1 ± 3.0
	HCT L/L	0.182 ± 0.018	0.333 ± 0.024	0.521 ± 0.030
	MCV fL	72.0 ± 5.0	80.5 ± 5.0	89.0 ± 5.0
	MCH pg	24.9 ± 2.8	29.0 ± 2.0	31.8 ± 2.0
	MCH fmol	1.54 ± 0.18	1.80 ± 0.16	1.97 ± 0.16
	MCHC g/dL	34.6 ± 3.6	36.0 ± 2.8	35.7 ± 2.8
	MCHC g/L	346 ± 36	360 ± 28	357 ± 28
	MCHC mmol/L	21.4 ± 2.3	22.3 ± 1.8	22.1 ± 1.8
	RDW-SD fL	36.0 ± 8.0	36.0 ± 8.0	42.0 ± 8.0
	RDW-CV %	13.5 ± 5.0	11.0 ± 5.0	11.0 ± 5.0
	PLT × 10 ³ /μL	62 ± 18	195 ± 30	392 ± 60
	MPV fL	9.2 ± 3.0	9.2 ± 3.0	9.5 ± 3.0
	W-SCR (LYMPH) %	12.0 ± 4.0	31.5 ± 4.0	62.0 ± 7.0
	W-MCR (MXD) %	19.5 ± 5.0	17.5 ± 4.0	12.5 ± 4.0
	W-LCR (NEUT) %	68.5 ± 7.0	51.0 ± 5.0	25.5 ± 5.0
	W-SCC (LYMPH) × 10 ³ /μL	0.3 ± 0.2	2.5 ± 0.6	12.2 ± 3.0
	W-MCC (MXD) × 10 ³ /μL	0.5 ± 0.3	1.4 ± 0.5	2.5 ± 1.1
	W-LCC (NEUT) × 10 ³ /μL	1.8 ± 0.5	4.0 ± 0.9	5.0 ± 1.7

Before using, refer to the instruction sheet for mixing directions.

For technical assistance in the USA and Canada call Technical Service at (800) 523-3395.

All brands and products are trademarks or registered trademarks of their respective companies.

R&D Systems, Inc.
614 McKinley Place NE
Minneapolis, MN 55413 USA
AS066-018 Rev. 11/17

INTENDED USE

CBC-ST Plus is a control designed to monitor values on automated and semi-automated impedance type hematology analyzers. It can also be used for manual methods. Please refer to the assay table for specific instrument models.

SUMMARY AND PRINCIPLE

It is an established laboratory practice to use a stable control to monitor the performance of diagnostic tests. This control is composed of stable materials that provide a means of monitoring the performance of hematology blood cell counters. It is sampled in the same manner as a patient specimen.

REAGENTS

CBC-ST Plus is an *in vitro* diagnostic reagent composed of human erythrocytes, simulated leukocytes, and mammalian platelets suspended in a plasma-like fluid with preservatives.




PRECAUTION

CBC-ST Plus is intended for *in vitro* diagnostic use only by trained personnel.



WARNING:

POTENTIAL BIOHAZARDOUS MATERIAL. For *in vitro* diagnostic use. Each human donor/unit used in the preparation of this product has been tested by a FDA licensed method/test and found to be negative or non-reactive for the presence of HBsAg, Anti-HCV, NAT testing for HIV-1, HCV (RNA) and HIV-1/2. Each unit is also negative by a serological test for Syphilis (RPR or STS). Because no test method can offer complete assurance that infectious agents are absent, this material should be handled as potentially infectious. When handling or disposing of vials follow precautions for patient specimens as specified in the OSHA Bloodborne Pathogen Rule (29 CFR Part 1910, 1030) or other equivalent biosafety procedures.

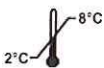

STABILITY AND STORAGE

Store CBC-ST Plus upright at 2 - 8° C (35 - 46° F) when not in use. **Protect vials from overheating and freezing.** Unopened vials are stable through the expiration date. Opened vials are stable for 14 days, provided they are handled properly.

INDICATIONS OF DETERIORATION

After mixing, product should be similar in appearance to fresh whole blood. In unopened tubes/vials, the supernatant may appear cloudy and reddish; this is normal and does not indicate deterioration. Other discoloration, very dark red supernatant or unacceptable results may indicate deterioration. **Do not use the product if deterioration is suspected.**



INSTRUCTIONS FOR USE

- Remove tubes/vials from the refrigerator and allow to warm to room temperature (15 to 30°C or 59 to 86°F) for 15 minutes before mixing.
- To mix, hold a tube/vial horizontally between the palms of the hands. **Do not pre-mix on a mechanical mixer.**
 - Roll the tube/vial back and forth for 20 - 30 seconds; occasionally invert the tube/vial. Mix vigorously, but do not shake.
 - Continue to mix in this manner until the red cells are completely suspended. Tubes/vials stored for a long time may require extra mixing.
 - Gently invert the tube/vial 8 - 10 times immediately before sampling.

- Analyze the sample as instructed in the Quality Control section of the Operator's Manual for your instrument.
- After sampling:
 - If tube/vial has been open for sampling, clean residual material from the cap and tube rim with a lint-free tissue. Replace the cap tightly.
 - Return tubes/vials to refrigerator within 30 minutes of use.

EXPECTED RESULTS

Verify that the lot number on the tube/vial matches the lot number on the table of assay values. Assay values are determined on well-maintained, properly calibrated instruments using the instrument manufacturer's recommended reagents. Reagent differences, maintenance, operating technique, and calibration may contribute to inter-laboratory variation.

PERFORMANCE CHARACTERISTICS

Assigned values are presented as a Mean and Range. The Mean is derived from replicate testing on instruments operated and maintained according to the manufacturer's instructions. The Range is an estimate of variation between laboratories and also takes into account inherent imprecision of the method and expected biological variability of the control material.

Assay values on a new lot of control should be confirmed before the new lot is put into routine use. Test the new lot when the instrument is in good working order and quality control results on the old lot are acceptable. The laboratory's recovered mean should be within the assay range.

For greater control sensitivity each laboratory should establish its own mean and acceptable range and periodically reevaluate the mean. The laboratory range may include values outside of the assay range. The user may establish assay values not listed on the Assay Sheet, if the control is suitable for the method.

LIMITATIONS

The performance of this product is assured only if it is properly stored and used as described in this insert. Incomplete mixing of a tube/vial prior to use invalidates both the sample withdrawn and any remaining material in the tube/vial.

TECHNICAL ASSISTANCE AND CUSTOMER SERVICE

For assistance in resolving control recovery problems, please call Technical Service at (800) 523-3395. For additional information on R&D Systems, Inc. hematology controls and calibrators, or to place an order, call Customer Service at (800) 428-4246.

QUALITY CONTROL PROGRAM

For information on CBC-Monitor, our Inter-Laboratory Quality Control Program, call (800) 523-3395 ext. 4435.

All brands and products are trademarks or registered trademarks of their respective companies.



R & D Systems, Inc.
614 McKinley Place NE
Minneapolis, MN USA 55413

IS066-007 Rev. 09/13



EUROCELL Diagnostics
19 Rue Louis Delourmel
35230 NOYAL CHATILLON/
SEICHE
France

